


SeDaCom

SEDACOM (Serial DATA COMMunication) is a Software that allows to display the data measured by different equipments LETICA, associated from communication RS 232 and provided without additional cost.

The gathered data are elaborated in ASCII format, being therefore applicable to most of statistical Software.

It fulfills GLP rules and increases the functionality of the equipment, saving and listing automatically all the data of the work sessions.

It can work simultaneously with until 10 equipments in the same PC, limited solely by the number of series ports available.

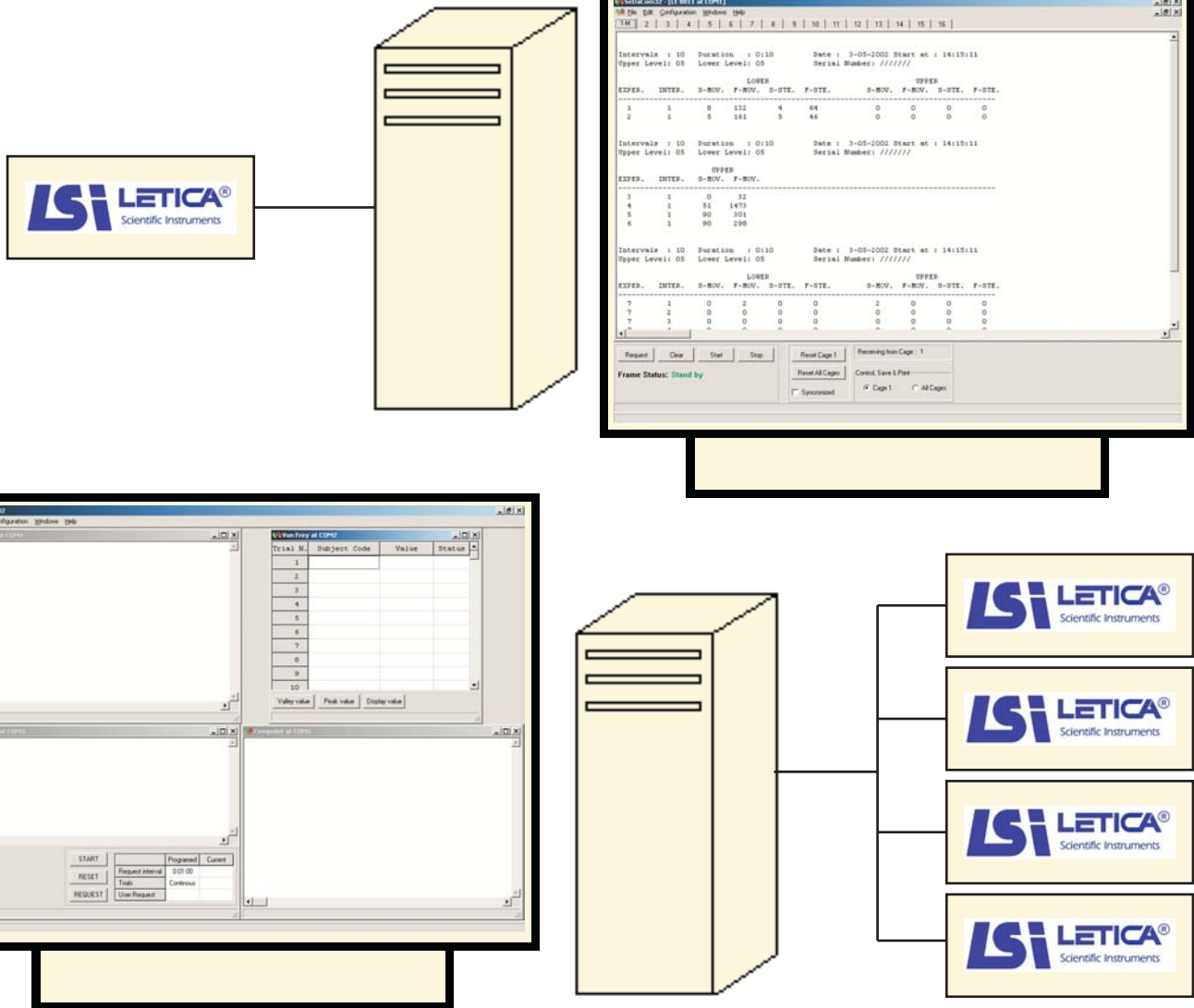

The LETICA's equipments that include SeDaCom are:

- *Plethysmometer*
- *Analgesymeter*
- *Non Invasive Blood Pressure Meter*
- *Rotometer*
- *Thermometer*
- *Grip Test*
- *Rota-Rod*
- *IR Actimeter*

EXPANSION CARDS (not included)

If the PC has two RS 232 ports, they can be connected directly two equipement. If it's necessary to connect more systems, they must be used Expansion Cards, that allow to increase the connection possibilities:

*Ref. SUNIX 405 5A for 4 ports more
SUNIX 406 5A for 8 ports more*

SeDaCom recognizes automatically the number of systems connected.

PC requirements: Each PC that runs at least with Windows 95.

SeDaCom is compatible with Windows 95, 98, ME, NT, 2000 and XP.